


HIGH-CAPACITY DEIONIZER WATER CARTRIDGE

Part Number – DB-10446


- The 'T bag' (Deionizer bag, P/N: DB-H200-031) is placed inside the filter cartridge instead of directly into the water tank. All the deionized water entering the electrolytic cell passes under pressure through this deionizer cartridge.
- **Efficient filtration of all the water** – this allows the generator operating parameters to be increased resulting in improved cell performance and reliability.
- **Easy to install and change** – with the quick connect fittings the cartridge is easy to install during the initial installation and simple to change when required.
- **Robust against challenging environments** – should the laboratory ambient air contain higher than normal levels of VOC's, the filter cartridge will help to maintain the water quality.

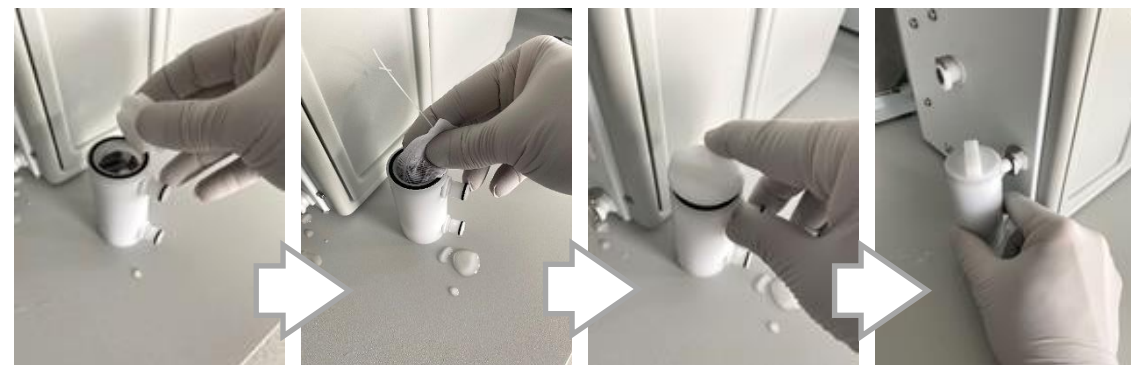

Deionizer bags

Part Number DB-H200-031

- Deionizer bags inside the foil packet. This protects the bags from ambient air, particularly CO₂, during storage and extends the shelf-life.
- There is 1 x bag in each foil packet
- Part number DB-H200-031 contains 2 x deionizer bags

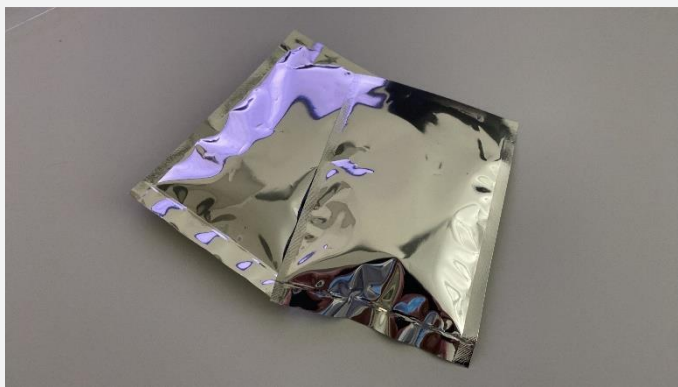

- Deionizer bags after removal from the foil packet.

